

### **Table of Contents**

HEPA-CARE® Patient Isolation Systems	3-5
HEPA-CARE® Replacement Filters & Accessories Chart	6
Monitors & Instruments	6
Surge Capacity Anteroom & Containment Modules	7
Isolation & Control Products For Construction & Renovation	8-9
PREDATOR® & HEPA-AIRE® Replacement Filters & Accessories Chart	9
Consumables & PPE	10-11
Contact Us	12

CANADA: 1-800-827-6443 • USA: 1-800-634-9091 • ABATEMENT.CA



## HEPA-CARE® Patient Isolation Systems

## HEPA-CARE® Patient Isolation Systems

For almost 30 years, HEPA-CARE® in-room HEPA filtration units have set the standard for reliable and economical solutions for infectious patient isolation. These unique systems are used today in thousands of facilities throughout North America. Two types of HEPA-CARE systems are available: portable models, which can be moved within the facility and deployed wherever they are needed; and, ceiling-mounted models, which mount into a drop-ceiling grid.

#### All models feature:

- Powder-Coated Aluminum Cabinet
   Professional, "hospital" appearance; easy to maintain.
- Strong, leak-resistant cabinet construction.
- Unique Zero Bypass Dual Cabinet Design Protection against contaminated air bypass.
- Intuitive capacitive touchscreen with 5 airflow settings.
- Password Protection
   Secure, tamper resistant controls.
- Energy Efficient EC motorized impellers Reduce energy use and noise.
- Anti-Microbial MERV-8 Pre-Filter
   Protects against microbial growth & degradation.
- Optional VAPOR-LOCK® Carbon Pre-Filter Removes low concentrations of gases and odours.
- All HEPA-CARE® HC500 series and HC800 series models are equipped with an H14 Type HEPA filter tested and certified to capture at least 99.995% of microscopic particles 0.3 microns in size



HC800FD PORTABLE AIR PURIFICATION SYSTEM

## NRTL Safety Certification

HEPA-CARE® systems are registered as Class II Medical Devices under FDA 510K, and certified by a Nationally Recognized Testing Laboratory (NRTL), to comply with relevant OSHA, UL and CSA safety standards.

## **Export Products**

Most of the Abatement Technologies electric-powered products in this brochure are available in 230 volt/50 hertz and 230 volt/60 hertz models for overseas customers. Contact our Healthcare Export Sales Department for more information.





## HEPA-CARE® Patient Isolation Systems

## HEPA-CARE® HEPA Filtration Systems

Thousands of health care facilities rely on HEPA-CARE systems from Abatement Technologies for infectious patient isolation. The two exceptionally mobile portable models can be quickly rolled to wherever needed in the facility and used to convert almost any room into a CDC-compliant patient isolation room. The ceiling models mount into a ceiling grid, completely out of the way of patients and staff. All four models are NRTL certified to comply with UL and CSA standards, and come equipped with: touch screen controls with password protection; five airflow settings; an anti-microbial treated pre-filter, and; a 99.995%-efficient HEPA filter.

#### HC800FD Portable Model



The HC800FD is the "premier standard" of portable systems, with top-of-the-line performance and features, and unmatched reliability. It can typically maintain 12+ ACH in rooms 4,000 to 4,625 cubic feet, depending on the installation.

- Peak rated airflow from 150cfm to 925cfm
- Louvered front air inlet & optional top & rear exhaust outlets
- Pharmaceutical grade HEPA with silicone-gel seal & anodized aluminum frame
- Energy Efficient EC motorized impellers reduce energy use and noise
- · Hinged, front door secures with two Phillips machine screws
- Clean and Sterilize the Air with optional Germicidal UV Disinfection. (HC800FD-UV)

#### PART NO: HC800FD

#### **HC500FD Portable Models**



This compact model is designed for areas that don't require as much airflow as the HC800FD. It can typically maintain 12+ ACH in rooms 2,250 to 2,375 cubic feet, depending on the installation.

- Peak rated airflow 150cfm to 475cfm
- Louvered front air inlet & optional top & rear exhaust outlets
- Pharmaceutical grade HEPA with silicone-gel seal & anodized aluminum frame
- Energy efficient EC motorized impeller
- Clean and Sterilize the Air with built-in optional Germicidal UV Disinfection (HC500FD-UV)

PART NO: HC500FD

#### HEPA-CARE® Germicidal UV Modules

Modules utilize a bank of four high-output UV lamps as a supplemental infection control measure.

#### UV800F UV Module



The UV800F is designed specifically for use with the HEPA-CARE® HC500FD and HC800FD portable units, mounting onto the top exhaust outlet. Air flowing through the UV800F is irradiated with a high dosage of germicidal UV (UV-C) energy.

- On/off switch so the portable unit can be run with or without UV-C
- Safety interlock switch protects against UV exposure

PART NO: UV800F



## HEPA-CARE® Patient Isolation Systems

## HC800CD Ceiling-Mounted Model



The HC800CD is similarly equipped to the HC800FD. It can typically maintain 12+ ACH in rooms as large as 3,625 to 4,250 cubic feet, depending on the installation.

- Peak rated airflow 725cfm 850cfm, depending on the installation
- Low-cost option to 'design & build' isolation rooms with dedicated HVAC systems
- Mounts into a 2 ft x 4 ft ceiling grid
- Pharmaceutical grade HEPA with silicone-gel seal & anodized aluminum frame
- Easy, in-room access for filter changes (no need to disturb ceiling tiles)
- Separate, wall mounted control box can be mounted outside the room or remotely

PART NO: HC800CD

## HC500CD Ceiling-Mounted Air Purification System



This compact model is designed for areas that don't require as much airflow as the HC800 series. It can typically maintain 12+ ACH in rooms 2,250 to 2,750 cubic feet, depending on the installation.

- Installs in a 2 ft. x 2 ft. drop ceiling grid
- Pharmaceutical grade HEPA with silicone-gel seal & anodized aluminum frame
- Easy, in-room access for filter changes (no need to disturb ceiling tiles)
- Separate, wall mounted control box can be mounted outside the room or remotely
- Five motor speed settings, with peak output of 200cfm up-to 550cfm\*
- Password protected to prevent tampering

Note: Use airflow ratings and air changes per hour (ACH) as a guideline only. Estimates are based on factory and independent testing @ 120 VAC with an air straightener and a traverse of readings taken with a computing vane-anemometer. Actual results may vary.

PART NO: HC500CD

## UV400C-PT & UV400C Ceiling-Mounted UV Modules





The UV400C-PT is designed specifically for irradiating exhaust air from any model. The UV400C is a module that can also be used in place of standard return-air grills in high infection risk areas.

- Mounts into a 2 ft. x 2 ft. drop ceiling grid
- Safety interlock shields protect against UV exposure

PART NO: UV400C-PT & UV400C

## HFPA-CARF® RFPI ACEMENT FILTERS & ACCESSORIES CHART

FILTER OR ACCESSORY	HC800FD	HC500FD	HC800CD	HC500CD	UV800F	UV400-PT	UV400C
Stage 1							
Anti-Microbial Pleated Prefilter	HC1802AM	HC1802AM	HC1802AM	HC1802AM	N/A	N/A	N/A
Optional Prefilter:							
VAPOR-LOCK® Carbon Filter	HC1802CF	HC1802CF	HC1802CF	HC1802CF	N/A	N/A	N/A
Final Stage							
99.97% HEPA Filter	HC1836-H14	HC1818-H14	HC1836-H14	HC1818-H14	N/A	N/A	N/A
UVGI Lamp	N/A	N/A	N/A	N/A	UV425	UV425	UV425
Clear Heavy-Duty	H2030-08	H2030-08	N/A	N/A	N/A	N/A	N/A
PVC Flex Duct, 8" x 25'							
Clear Wire Reinforced	H2030-8-12	H2030-8-12	N/A	N/A	N/A	N/A	N/A
PVC Flex Duct, 8" x 12'							
Corrugated Aluminum	H20508-8-1	H20508-8-1	H20508-8-1	N/A	H20508-8-1	H20508-8-1	H20508-8-1
Semi-Rigid Duct, 8" x 8'							
Quick-Connect Collar	H2240-08K	H2240-08K	H2240-08K	H2240-08K	N/A	N/A	N/A
Flex Duct Attachment Clamp	H2350-08						



## Monitors and Instruments

## Single Room Pressure Monitor



PART NO: RPM-RT1 OR RPM-RT2

Abatement Technologies® RPM Series Differential Pressure Monitors provide 24/7 monitoring and recording of differential air pressure relationships between multiple rooms or areas within a healthcare facility for patient isolation environments.

- User-friendly, colour touch screen display
- Differential pressure measurement, ranges of -0.10inWC (-25 Pa) to +0.10inWC (+25 Pa)
- Archived logs are downloadable in CSV format to a USB flash drive
- Password protection feature (can be disabled)
- Visual & audible alarms with audio alarm delay & timed shutoff features
- (Compatible with) BACnet and N2 building automation systems
- Optional Temperature and Relative Humidity (RPM-SENSOR)

## Portable Differential Pressure Monitor



PART NO: PPM3-S

The PPM3-S Portable Room Pressure Monitor from Abatement Technologies® is an essential tool for measuring the realtime pressure differential between two adjacent areas within a facility.

- User-friendly, colour touch screen display
- Differential pressure measurement, ranges of -0.10inWC (-25 Pa) to +0.10inWC (+25 Pa)
- Archived logs are downloadable in CSV format to a USB flash drive
- Password protection feature (can be disabled)
- Visual & audible alarms with audio alarm delay & timed shutoff
- 6VDC @ 1 Amp output power supply/ charger with 10' cable
- Up to 6 months of battery life operating in battery mode.
- Protective carrying case



## Surge Capacity Anteroom & Containment Modules

## Temporary, Reusable, Adaptable Modular Containment Kit



The AGSHIELD-AR1 starter package includes a 4ft outward swinging door, a 4ft panel, a 2ft panel, a 2ft panel with an exhaust port, 2, 90° Fixed Corners, 2 Wall Seal Drywall adapter kits, closed cell foam gasket, all necessary hardware and a transportation Panel Cart. Over 40 configurations are possible using this kit.

- Easily create different size anterooms or other configurations with one simple kit
- ICRA compliance up to Class 4
- Quick assembly with a single tool
- Secure environment with locking doors and exhaust plate
- Clean, professional, high-end appearance

PART NO: AGSHIELD-AR1

## ISOOUAD® 10810 Isolation Tent



**PART NO: 10810** 

The IQ810 is a flexible, inexpensive and easy way to quickly increase surge patient isolation capacity that one person can pop up in seconds.

- Fits almost any bed or gurney
- Use with a HEPA filtration device to create negative pressure within the module
- Pulls in "clean air at the foot of the bed & exhausts HEPAfiltered air behind the patient's head

## ANTE UP® HC7917 Inflatable Emergency Anteroom



PART NO: HC7917

The HC7917 can be mounted inside a patient room doorway or outside the door, entirely in the corridor.

- Inflates in minutes with the included 110-volt compressor/inflator
- Fits any 46" 48" wide, 80" tall doorway & can be expanded for double doors
- Large zippered doors at either end are wide enough for a gurney
- Large windows with roll-up covers at each end
- Portable air scrubber & HEPA vacuum connection ports
- Tough fabric meets relevant fire/life safety standards

## AIRE GUARDIAN® AG3100TMK Temporary Transition/Anteroom Module



PART NO: AG3100TMK

The AG3100TMK includes an aluminum frame and a PVC enclosure. An optional transport cart can be purchased so the module can be quickly wheeled to a patient room, removed from the cart and sealed around the door frame.

- Enclosure flange fits any 46" 48" wide, 80" tall doorway
- Large zippered doors at either end are wide enough for a gurney
- Large view windows
- Designed for use with a HEPA-CARE in-room HEPA filtration unit
- Tough military tent fabric meets relevant fire/life safety standards



# ISOLATION & CONTROL PRODUCTS FOR CONSTRUCTION & RENOVATION

### PREDATOR® Portable Air Scrubbers

Abatement Technologies® developed the original portable air scrubber more than 20 years ago and our products continue to set the standard for quality, performance, portability, reliability and ease of use.

## PREDATOR® 750 Portable Air Scrubber



- Airflow: Variable speed, rated 200cfm -750cfm
- Six (6) air changes per hour (ACH) for 5,000 - 6,000 cubic foot spaces
- Weight: 35 pounds with filters
- Cabinet: Durable, Polyethylene Made From UL94HB Flame Retardant Resin
- Peak power draw: 1.9 amps or less, daisychain up to 5 units per 15A circuit

PART NO: PRED750

## PREDATOR® 1200 Portable Air Scrubber



PART NO: PRED1200

- Airflow: Variable speed, rated 300cfm -900cfm
- Six (6) air changes per hour (ACH) for 6,000
  9,000 cubic foot spaces
- Weight: 70 pounds with filters
- Built-in, two-wheel dolly and durable handle for maximum portability
- Cabinet: Durable, Polyethylene Made From UL94HB Flame Retardant Resin
- 3 units can be daisy-chained on one 15-amp circuit

## HEPA-AIRE ® Portable Air Scrubbers

HEPA-AIRE® models feature tough, corrosion-resistant metal cabinets fabricated using airtight "zero bypass" rivet construction.

#### HEPA-AIRE® 750 Portable Air Scrubber



PART NO: PAS750

- Airflow: Variable speed, rated 200cfm - 750cfm
- Six (6) air changes per hour (ACH) for 6,000 – 7,500 cubic foot spaces
- Weight: 33.5 pounds with filters
- Cabinet: Anodized aluminum
- Peak power draw: 1.9 amps or less, daisy-chain up to 5 units per 15A circuit

## HEPA-AIRE® 1200 Portable Air Scrubber



PART NO: PAS1200

- Airflow: Variable speed, rated 300cfm - 900cfm
- Six (6) air changes per hour (ACH) for 7,500 - 9,000 cubic foot spaces
- Weight: 85 pounds with filters
- Upright model, with nonmarking 8" rear wheels & front swivel casters
- Cabinet: Stainless steel
- 3 units can be daisy-chained on one 15-amp circuit

## HEPA-AIRE® 1600 High-Static Portable Air Scrubber



- Airflow: One speed, rated 1,900cfm
- Six (6) air changes per hour (ACH) for 17,500 - 19,000 cubic foot spaces
- Weight: 131 pounds
- Cabinet: Stainless steel
- Peak power draw: 8.5 amps or less

PART NO: PAS1600SHS



## ISOLATION & CONTROL PRODUCTS FOR **CONSTRUCTION & RENOVATION**

## HEPA-AIRE® 2400 Ultra-High Performance Portable Air Scrubber



- scrubber available
- Six (6) Air changes per hour (ACH) for 18,000 21,000 cubic foot spaces
- Weight: 170 pounds
- Cabinet: Stainless steel, with detachable top module
- Peak power draw: 12 amps or less
- Three standard filtration stages: 1" MERV 4 natural fiber filter, 2" MERV 8 pleated filter & HEPA filter

• Airflow: Two speed, rated 700 cfm (low)/ 2,100 cfm (high), the highest performance 120V/15A air

- Optional pre-filter: 2" 500gsm VAPOR-LOCK® carbon filter (adds 4th filtration stage)
- Each pre-filter has its own shelf to minimize bypass

PART NO: PAS2400

## HEPA-AIRE® 5000 High Performance Portable Air Scrubber



PART NO: PAS5000

- Airflow: One speed, rated 4,000 cfm, the highest performance air scrubber available
- Six (6) air changes per hour (ACH) for 35,000 40,000 cubic foot spaces
- Weight: 332 pounds
- Cabinet: Galvanized 20 gauge steel
- Peak power draw: 23 amps or less on a 230V/30A/60Hz electric supply
- Two standard filtration stages (2 of each filter required): 2" MERV 8 pleated pre-filters & **HEPA** filters
- Optional pre-filters (2 required): 2", 500gsm VAPOR-LOCK® carbon filters (replace standard pre-filters)
- Pre-filters mount in door for protection against bypass

## PREDATOR® & HEPA-AIRE® Replacement Filters & Accessories Chart

FILTER OR ACCESSORY	PRED750	PRED1200	PAS750	PAS1200	PAS1600SHS	PAS2400	PAS5000
Stage 1: Large Particle Filter	F621	F621	F621	F621	H2003	F1821	H2001(2)
Stage 2:							
MERV 8 Pleated Particle Filter	H502	H502	H502	H502	H2002	H1902	H2002(2)
VAPOR-LOCK® Carbon	VL1002	VL1002	VL1002	VL1002	VL2024	VL2002	VL2024(2)
Stage 3: VAPOR-LOCK® Carbon	N/A	VL1002	N/A	VL1002	N/A	VL2002	N/A
Final Stage: 99.97% HEPA Filter	H161606-99	H610C-99	H161606-99	H610C-99	H24240699DG	H2418-99	H242406-99(2)
Lay Flat Outlet Ducting, 500' Roll	H2015-12	H2015-12	H2015-12	H2015-12	H2015-12	H2015-12	H2015-12
Clear Wire-Reinforced Mylar, 25'	H2020-10	H2020-10	H2020-10	H2020-10	H2020-12	H2020-12	H2020-12
Hvy Duty Wire-Reinforced PVC	H2030-10	H2030-10	H2030-10	H2030-10	H2030-12	H2030-12	H2030-12
Flex Duct Attachment Clamp (intake):	H2350-10	H2350-10	H2350-10	H2350-10	H2350-12	H2350-12	H2350-12(2)
Flex Duct Attachment Clamp (Exhaust):	H2350-10	H2350-10	H2350-10	H2350-10	H2350-12	H2350-12	H2350-14(2)

Note: Only bold-faced items come with the base unit. All others must be purchased separately.



## Consumables and PPE

In addition to the trusted products we make and sell, we have partnered with some of Canada's best brands to provide you with a more complete line of products to try and cover everything you do on site. Making Abatement Technologies® your single destination for quality equipment and supplies.

## PERSONAL PROTECTIVE EQUIPMENT



HAND PROTECTION:

Nitrile, Rubber, Cloth and Cut Resistant Work Gloves



**HEAD PROTECTION:** 

Certified Hard Hats and Accessories



**RESPIRATORS AND FILTERS:** 

Full Facepiece, Half Facepiece, Particulate and Combo Filters



DISPOSABLE DUST MASKS



EYE WEAR:

Safety Glasses and Chemical Splash Goggles



**FALL PROTECTION:** 

Harnesses, Lanyards, Self-Retracting Lifelines and Anchors



FIRST AID PRODUCTS:

Eyewash Stations, Complete First Aid Kits



DISPOSABLE PROTECTIVE CLOTHING:

Breathable, Non-Breathable Coveralls, Boot and Shoe Covers

## TOOLS AND ACCESSORIES



GRILL MASK AND DISPENSERS



STICKY FLOOR MATS



AIR SCRUBBER FILTERS



**TAPES** 



**DUCTING AND CLAMPS** 



Partner Brands

























## CONTAINMENT PRODUCTS



POLY SHEETING, CLEAR, RIP-PROOF, FIRE RETARDANT



FLOOR PROTECTION: Hard Floor, Carpet and Reusable Floor Runners



Disposable Zippers, Spring Loaded Poles, Door Kits



POLY HANGERS AND GRABBERS



Asbestos and Garbage. Available in Multiple Sizes and Thicknesses

**DISPOSABLE BAGS:** 



DISINFECTANTS




**DEGREASERS** 

## **ADHESIVES**



SPRAY GLUES

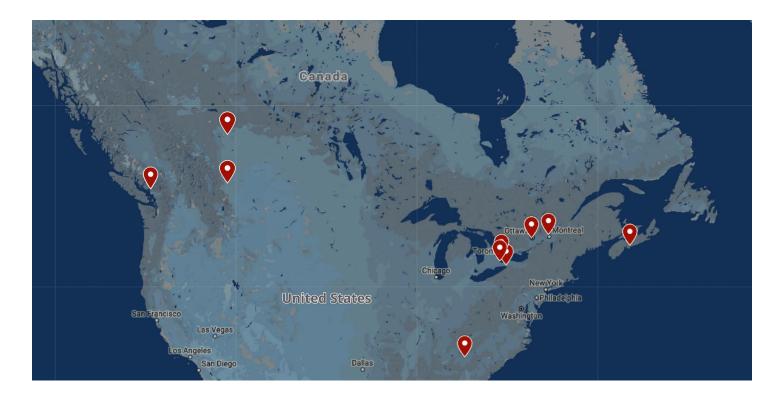

## Main Offices

Canada

Abatement Technologies Limited

7 High Street

Fort Erie, ON L2A 3P6

CAN. 800-827-6443

Intl. +01-905-871-4720

iaqinfo@abatement.ca

US

Abatement Technologies, Inc. 605 Satellite Blvd. Suite 300

Suwanee, GA 30024

U.S. 800-634-9091

Intl. +01-678-889-4200

iaqinfo@abatement.com

## **Branch Offices**

Vancouver

F102-19720-94A Avenue, Langley, BC V1M 3B7

778-372-2611

Edmonton

16403 - 117 Ave.

NW Edmonton, AB T5M 3V3

780-665-5165

Toronto

144 Sharer Rd. Vaughan, ON L4L 8P4

905-266-1322

Montreal

1950 Bernard Lefebvre Laval, QC H7C OA5

514-375-0606

Calgary

F-7004 5th St. SE Calgary, AB T2H 2G3

403-351-2223

Hamilton

428 Millen Rd. Unit 10 Stoney Creek, ON L8E 3N9

905-664-6075

Ottawa

51 Capital Drive Ottawa, ON K2G 0E7

613-749-9667

Dartmouth

121 IIsley Ave Unit 12 Dartmouth, NS B3B 1J1

902-334-0900

1.800.827.6443

ABATEMENT.CA IAQINFO@ABATEMENT.CA

© 2020, Abatement Technologies, All Rights Reserved. No reproduction of any portion of this advertisement is permitted without the express written permission of Abatement Technologies, All illustrations and specificatons contained in this advertisement are based on the latest product information available at the time of printing, Abatement Technologies reserves the right to (1) make changes at any time, without notice, to product designs, specifications, materials, and colors, (2) to introduce new models to the marketplace and, (3) discontinue products as it sees fit.